

**NOTES:**

- ① GAS OUTLETS TO BE MSC STANDARD BRAND OUTLET WITH CONNECTION TO MATCH EXISTING HOSPITAL EQUIPMENT.
- ② SERVICES REVERSED ON OPPOSITE HAND UNITS.
- ③ CEILING HEIGHT: 7'6" MIN TO 9'6" MAX. OTHER HEIGHTS ARE SUBJECT TO AN ADDITIONAL CHARGE. CONTACT YOUR SALES REP FOR DETAILS.
- ④ TWO SINGLE SIDED UNITS REQUIRED FOR BACK TO BACK INSTALLATIONS. TOTAL WALL THICKNESS WITH SINGLE LAYER OF 5/8" GYPSUM = 7.25".
- ⑤ APPROXIMATE JUNCTION BOX LOCATIONS.
- ⑥ MEDICAL GAS CONNECTIONS @ 85". EXACT HORIZONTAL SPACING TO BE DETERMINED IN SUBMITTALS.
- ⑦ INTERIOR STUD FRAMING 2-3/4" STUD WITH PRE-ATTACHED Z-BRACKET TO FIT A 3-5/8" METAL STUD WALL WITH 5/8" GYPSUM BOARD FINISH.
- ⑧ NO DRYWALL FASTENERS ARE TO BE USED IN THIS AREA.
- ⑨ PANEL LAMINATE COLOR TO BE DETERMINED. POWDER COATED SURFACES ARE TO MATCH FORMICA #923 "SURF".



**MODULAR SERVICES COMPANY** 500 EAST BRITTON ROAD  
OKLAHOMA CITY, OK 73114

REVISED

PROJECT

DUAL SILHOUETTE®  
SEMI-RECESSED  
MODEL 8600 (SP)


SCALE: 1" = 1'

DATE \_\_\_\_\_ DRAWN BY \_\_\_\_\_ DRAWING #  
PB-196

**LEGEND**

REQ	SYM	DESCRIPTION
3	□	OXYGEN OUTLET
2	A	AIR OUTLET
4	V	VACUUM OUTLET
4	E2	DUPLEX RECEPTACLE, 8300 IVORY, 20A, 125V, H.G., NORMAL POWER
6	E4	DUPLEX RECEPTACLE, 8300 RED, 20A, 125V, H.G., EMER. POWER
1	NC	NURSE CALL PROVISION
1	CB	CODE BLUE PROVISION
1	MO	MONITOR OUTPUT PROVISION
1	MC	MONITOR CHANNEL ONLY
8	HC	HORIZONTAL ACCESSORY CHANNEL

QTY	ACCESSORIES NOT SHOWN (PER UNIT)

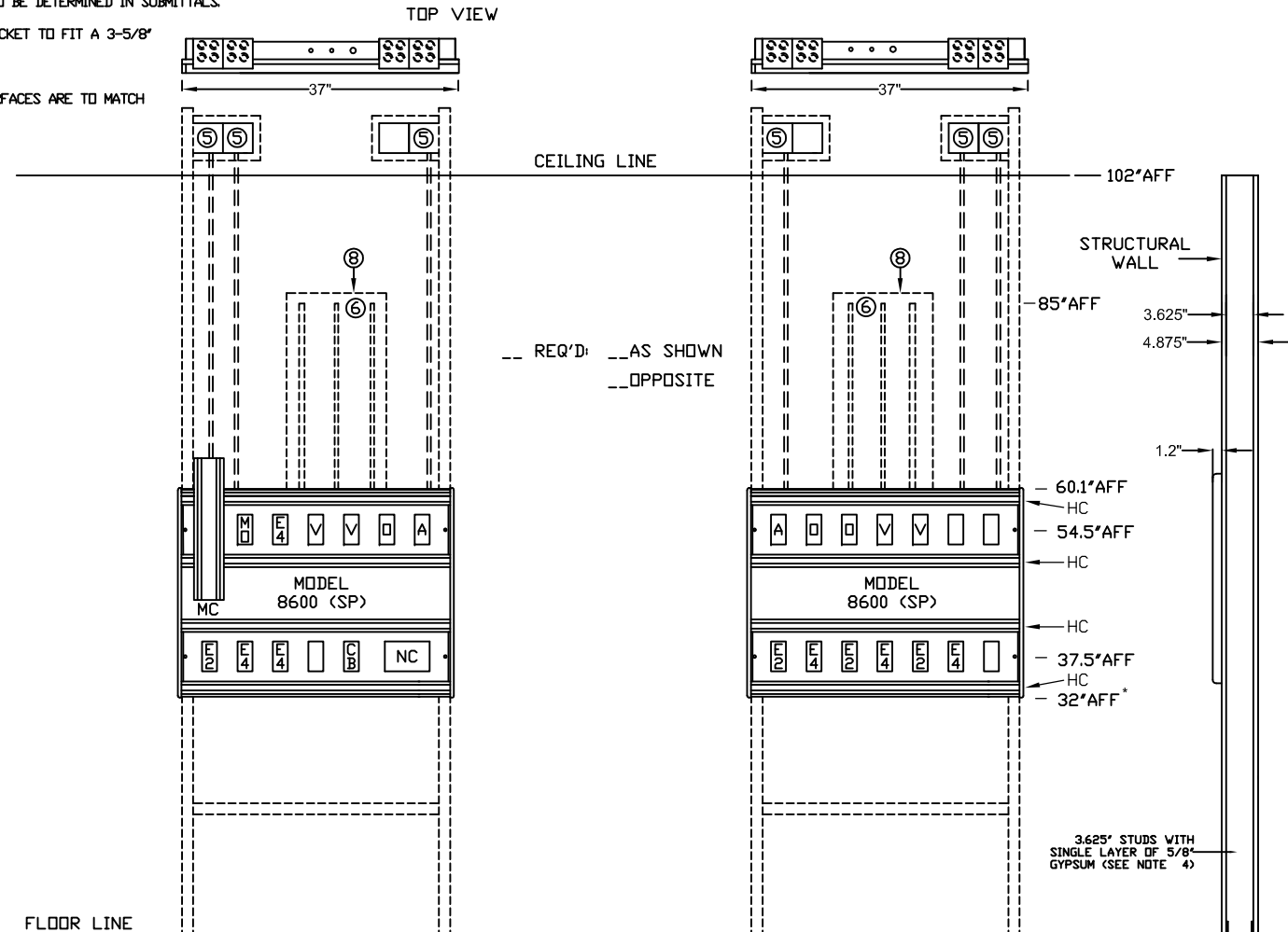

FLOOR LINE  
\* AFF (ABOVE FINISHED FLOOR)