

NOTES:

- ① JUNCTION BOXES: J1--NORMAL POWER, J2--EMERGENCY POWER, J3--NURSE CALL/CODE BLUE/TELEPHONE/MULTI PIN/MONITOR OUTPUT.
- ② MEDICAL GAS PIPING: OXYGEN & AIR--1/2" I.D., VACUUM 3/4" I.D.
- ③ INTERIOR STUD FRAMING IS TO FIT A 3-5/8" METAL STUD WALL WITH ONE LAYER OF 5/8" GYPSUM BOARD FINISH.
- ④ PANEL LAMINATE COLOR: TO BE DETERMINED. POWDER COATED & KYDEX SURFACES ARE TO MATCH 'SURF' OR 'PLATINUM GREY'.
- ⑤ EMERGENCY POWER, NORMAL POWER, AND LOW VOLTAGE FIELD CONNECTIONS ARE REQUIRED BETWEEN SECTIONS A&B AND B&C.

**MODULAR SERVICES COMPANY** 500 EAST BRITTON ROAD  
OKLAHOMA CITY, OK. 73114

PROJECT


REVISED

SILHOUETTE®  
SEMI-RECESSED  
MODEL 8600B

SCALE: = 1'


DATE \_\_\_\_\_ DRAWN BY \_\_\_\_\_ DRAWING #  
PB-169

SEE DRAWING LEG169 FOR WIRING DIAGRAM AND NOTES.

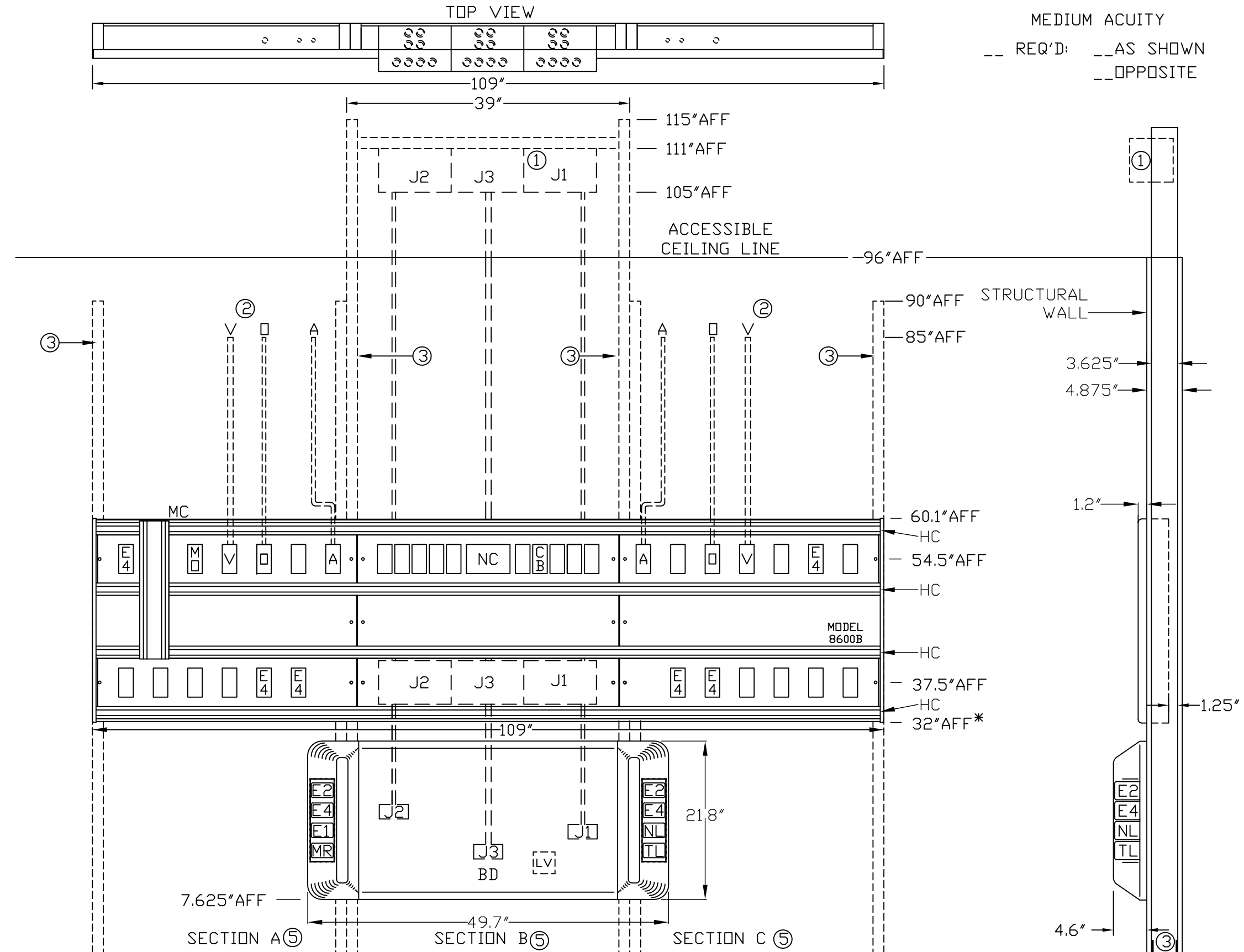

MEDIUM ACUITY  
-- REQ'D: -- AS SHOWN  
-- OPPOSITE

\* AFF (ABOVE FINISHED FLOOR)

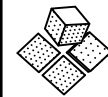


PROJECT

SCALE: =1'

DATE \_\_\_\_\_ DRAWN BY \_\_\_\_\_ DRAWING #  
LEG169

**LEGEND**

REQ	SYM	DESCRIPTION
2	□	OXYGEN OUTLET
2	A	AIR OUTLET
2	V	VACUUM OUTLET
1	E1	SINGLE RECEPTACLE, 8310 IVORY, 20A, 125V, H.G., NORMAL POWER
8	E4	DUPLEX RECEPTACLE, 8300 RED, 20A, 125V, H.G., EMER. POWER
2	E2	DUPLEX RECEPTACLE, 8300 IVORY, 20A, 125V, H.G., NORMAL POWER
1	NC	NURSE CALL PROVISION
1	CB	CODE BLUE PROVISION
1	MO	MONITOR OUTPUT PROVISION
4	HC	HORIZONTAL ACCESSORY CHANNEL
1	BD	BED DOCKER
1	TL	TELEPHONE OUTLET SINGLE CAT 5, RJ-45 RECEPTACLE
1	MR	MULTI-PIN BED RECEPTACLE 37 PIN CONNECTOR W/CABLE
1	NL	NIGHT LIGHT, SIERRA 2156, 4 WATT, 115V, LOUVERED FACE
1	MC	518-2302 MONITOR CHANNEL ONLY, PROVIDER™ CHANNEL MOUNT
1	LV	LOW VOLTAGE SWITCHING PACKAGE

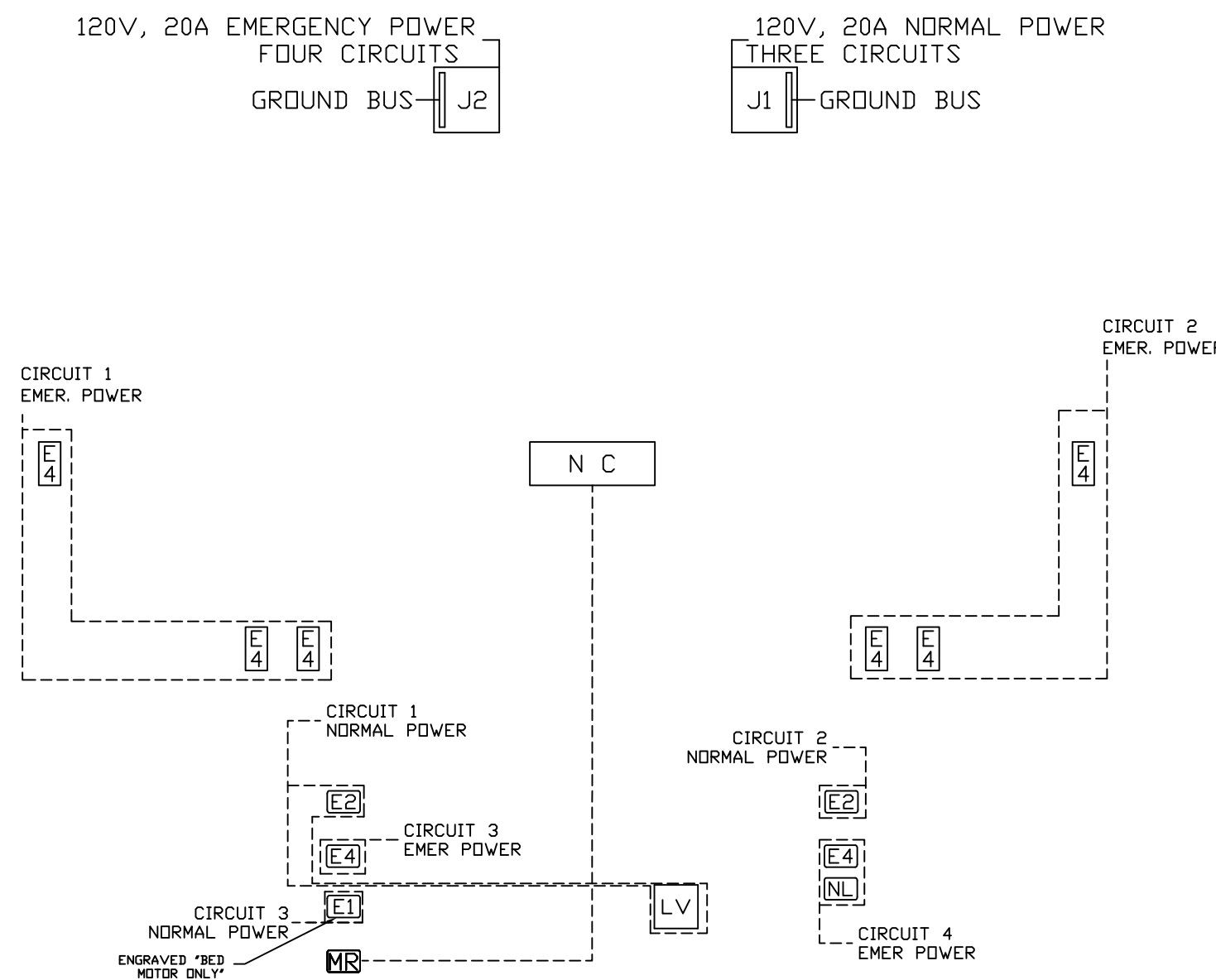

**WIRING DIAGRAM**

GROUND WIRING TO BE #10  
GREEN STRANDED COPPER.