




MODULAR
SERVICES COMPANY®

Quick Ship Solutions

MOBILE HEADWALL

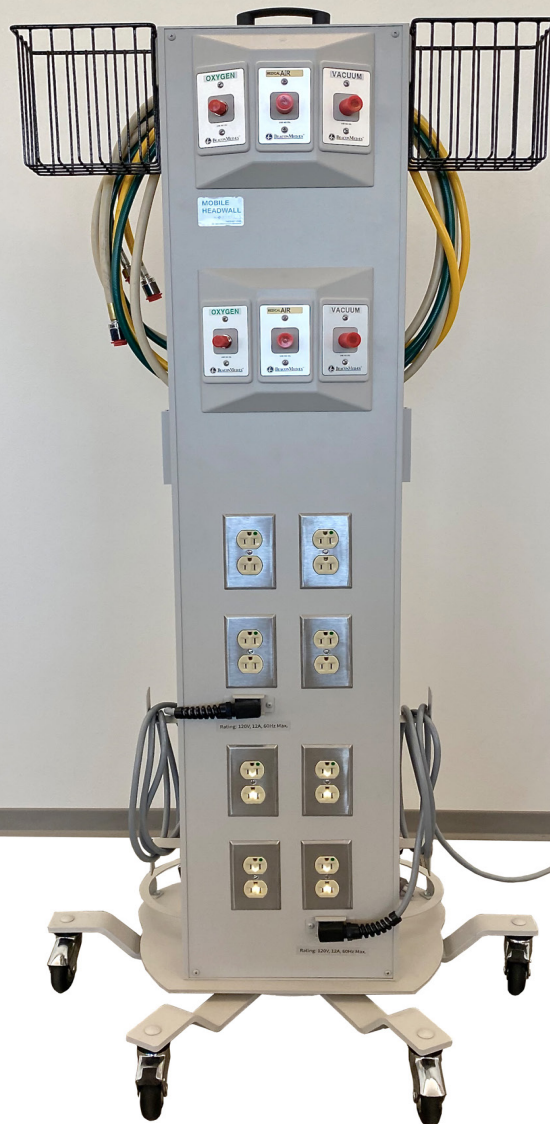


Designed for Emergency Response.

The Mobile Headwall provides a plug and play solution for medical gas and electrical services in temporary patient treatment areas. With a convenient top grab handle, and a mobile base, each unit can be quickly deployed to healthcare facilities for surge capacity needs.

MOBILE HEADWALL





Pre-mounted Vacuum Bottle Slides

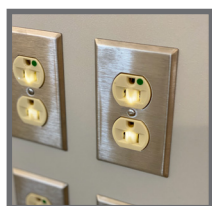

LED Indicator lights ensure power is obtained when plugged in to receptacle



Built-in cord management

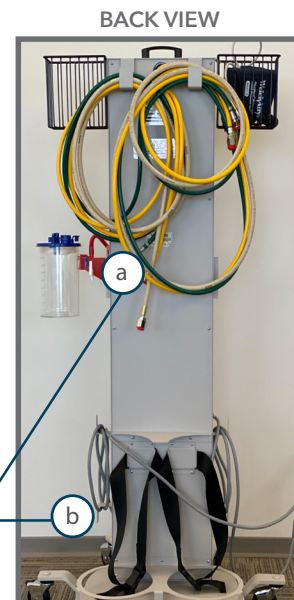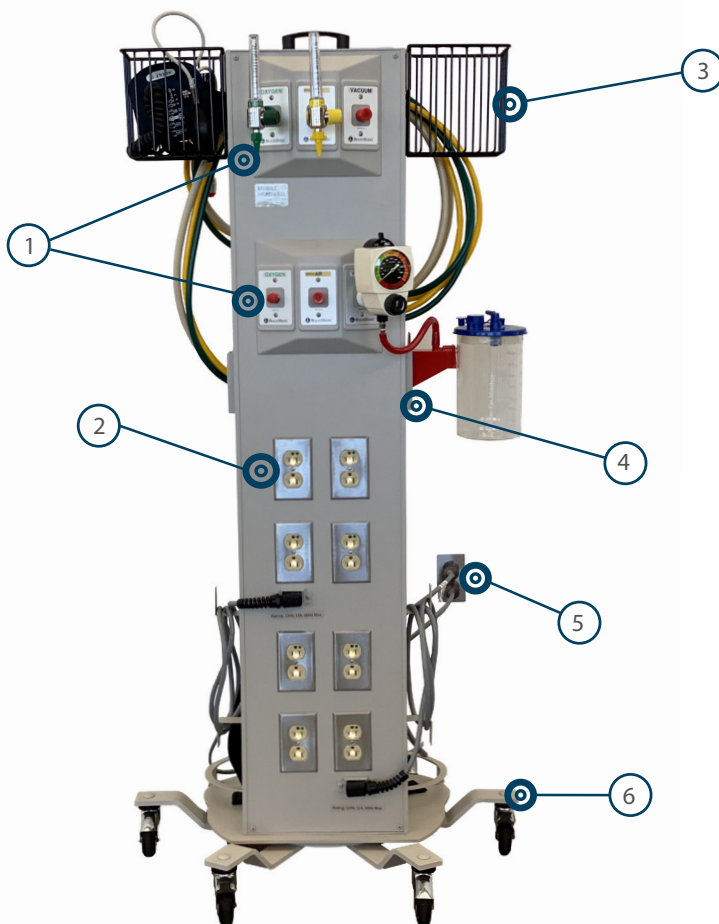

Product Features

- ① Low Acuity version (Model 10110) includes 3 Medical Gas outlets at 54" AFF, High Acuity version (Model 10120) includes 6 Medical Gas outlets, one tier at 42" AFF and one at 54" AFF
- ② Up to 2 Electrical Circuits: Low Acuity version includes 2 electrical outlets at 24" AFF, and 2 at 30" AFF. High Acuity version includes an additional 4 outlets, 2 at 10" AFF and 2 at 16" AFF
- ③ Baskets for Cuff and Supplies: 1 Basket included on Low Acuity version, 2 on High Acuity version
- ④ Pre-mounted Vacuum Bottle Slides: 1 on Low Acuity version, 2 on High Acuity version
- ⑤ Hospital grade, 10' power cord plugs into any standard 120-volt hospital grade power receptacle, normal or emergency
- ⑥ Mobile Base with 6 casters (3 locking) for added stability
- ⑦a Convert existing private rooms to multi-capacity by connecting the gas hoses at the back of the Mobile Headwall to existing medical gas outlets
- ⑦b Convert temporary spaces into patient treatment areas by securing 2 compressed gas cylinders to the back of the Mobile Headwall via cylinder housing and 2 straps